

# Burette Set

## 150mL Infusion Line

### DESCRIPTION

The Burette with Needleless Luer Access has a 150mL capacity vented drip chamber at 60 drops/mL, with a floating safety shut off valve. There is a vented spike, roller regulator and a swabable needleless access site at burette. While the tubing has an additional NeutralSite™ Neutral displacement valve, roller regulator, click clamp and a male luer lock connector with rotating collar.

### FEATURES & BENEFITS

- Vented perforator
- 2 x Rigid Roller Regulator for accurate flow
- NeutralSite™
- 150mL Burette
- Floating safety shut off valve
- Vent and needleless access site on burette chamber
- Tubing with click clamp
- Male luer lock with rotating collar
- Peel open pouch
- Single use only


### PRODUCT DATA

- Burette capacity: 150mL
- Tubing length: 195cm
- Total length: 230cm
- Tubing diameter: 3mm
- Drip chamber: 60 drops /mL
- Priming volume: 21.0mL
- NeutralSite™: Neutral displacement swabable needleless access site
- Sterile EO
- ARTG Number: 304860
- The product is MRI and CT suitable
- Not to be used with blood or blood products
- Not made with DEHP
- Not made with Natural Rubber Latex
- Material: DEHP Free PVC tubing, PETG, PP & Silicone, POM
- NeutralSite™ Materials: PP, ABS and Silicone
- Shelf Life: 3 Years

### ISO STANDARDS

- ISO 8536-5-2004 Infusion Equipment for Medical Use- Part 5: Burette infusion sets for single use, gravity feed
- ISO 8536-4-2010 Infusion Equipment for medical use - Part 4: infusion sets for single use, gravity feed

### PRODUCT DIAGRAM

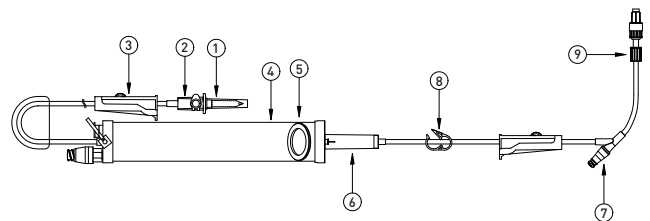

ITEM	DESCRIPTION
1	Vented Spike
2	Air Filter Cap
3	Roller Regulator
4	Burette 150mL
5	Stop Afloat
6	Chamber
7	NeutralSite™
8	Clamp
9	Male Luer Lock

### Burette Set

STERILE

RE-ORDER	DESCRIPTION	LENGTH	CARTON/QTY
IV040002	Burette Set 150mL 60drops/mL Infusion Line	230cm	25